

Distinguished Safety

Large Full-Flat working surface for optimal and safe handling of things placed in the chamber

Flexibility

Standard microscope provides more flexibility and ergonomic use to the SmartStation

Eyepoint Adjuster and Long Tilting Trinocular provides unrivalled ergonomics



The height of working surface is adjustable from 840 ~ 1090 mm



Mobility

Ergonomic handle bars on each side for comfortable mobility

Optional bracket for CO₂ gas tank is available

CO ₂ range:	0-10%
CO ₂ accuracy	+/-0.2%
Temp accuracy:	+/-0.5C
Temp uniformity:	+/-1.0C

Ex Dim (mm):	840 x 620 x 1190~1440
Weight:	80kg



Please contact us below
Email : info@astec-bio.com
TEL : +81-92-935-5666
Or your local ASTEC distributor

www.astec-bio.com



SmartStation for Smart Lab

Gives the Best to both laboratory and embryologists


ADVANCED TECHNOLOGY AT THE CORE

A WorkStation with Incubator-like Stable Environment helps to minimize the Stress on precious cells right from Retrieval room

Grade A Air with integrated HEPA and optional VOC Filters

Cell Manipulation or Observation in CO₂-controlled Optimal Environment


Integrated Flush-fitted, Heated Glass with Full-Flat working surface

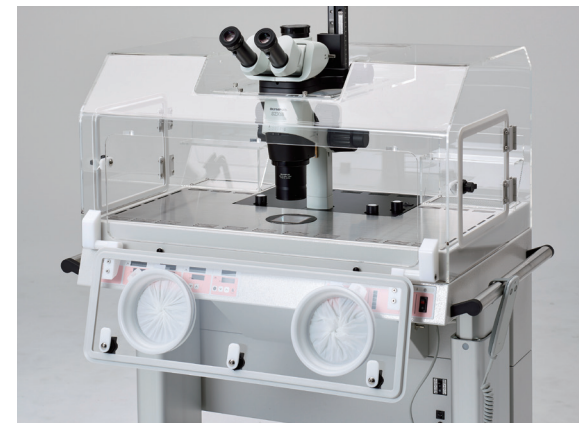

Front-Panel folds down for Easy and Safe Handling of Loaded ET Catheter



Comfortable, Full-Flat working area with embedded glass warmer and microscope illumination control can be fitted with most of the stereo microscopes available in the market

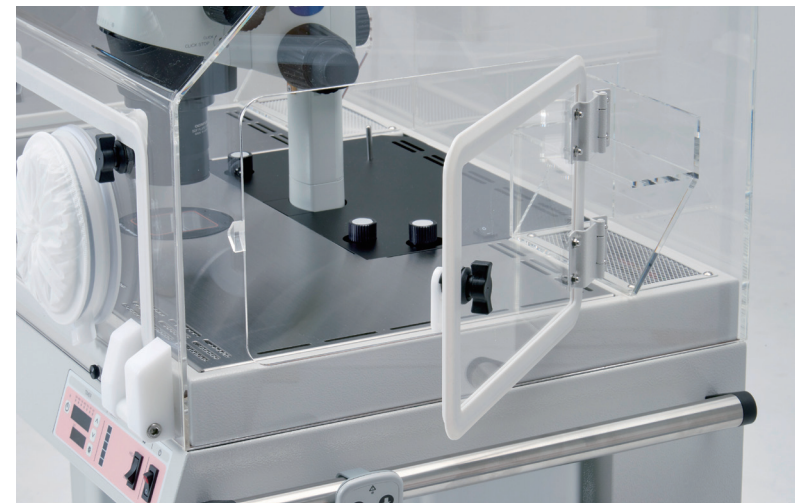


Large Side-doors on both sides for comfortable In-Out of Materials




Ergonomic Usability

A motorized stand helps you choose the most comfortable height to work with

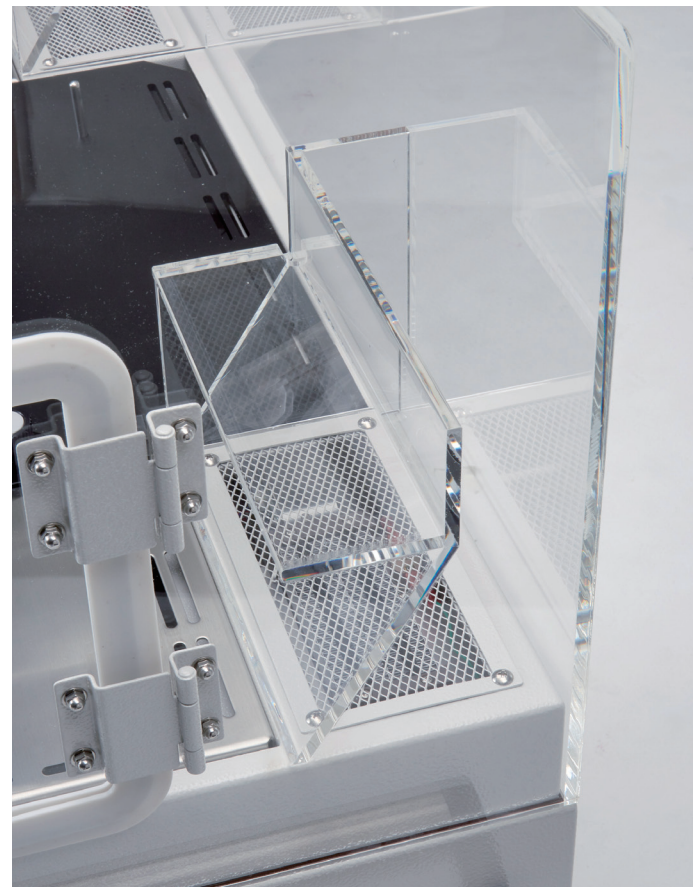

Utility SmartShelf on each side for those space-requirement without compromising with Table-top Working-Space

Control of Light and Vision Fields is all handy inside the working-chamber
Maintaining Hygiene at the highest level

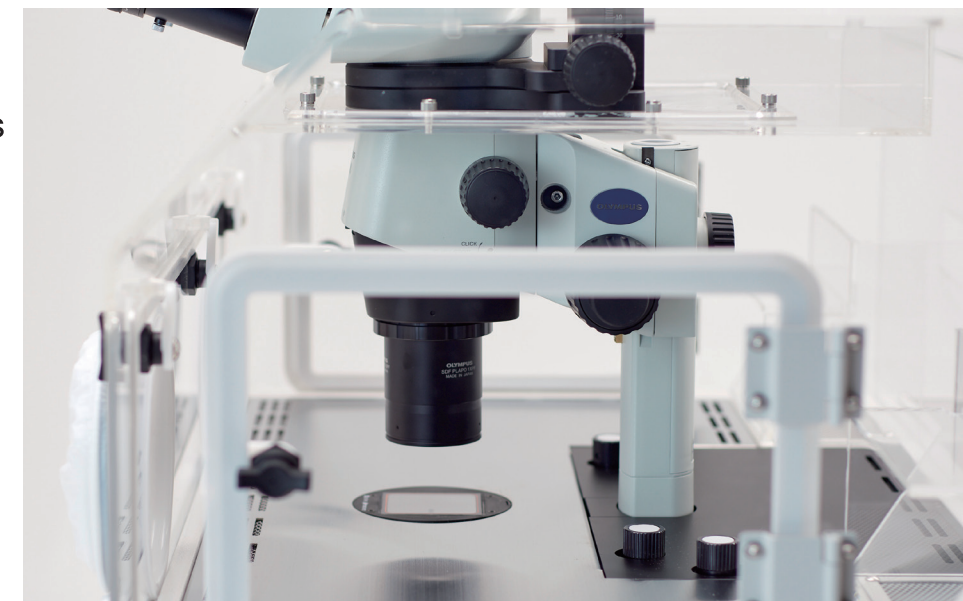

Designed to work with you