

ASTECC Gas line Filters


ASTECC-6: Six-month Filter

ASTECC-4: Four-month Filter



Achieve Better results with Gas Purified by Astec

Astec Gas Line Filters use a combination of activated carbon granules, a HEPA filter and pre-filters to remove VOCs and trap microorganisms, aldehydes and particles. Unique Micro Carbon Technology provides a large activated carbon surface area, enabling uniform gas flow and extremely high purifying capacity.

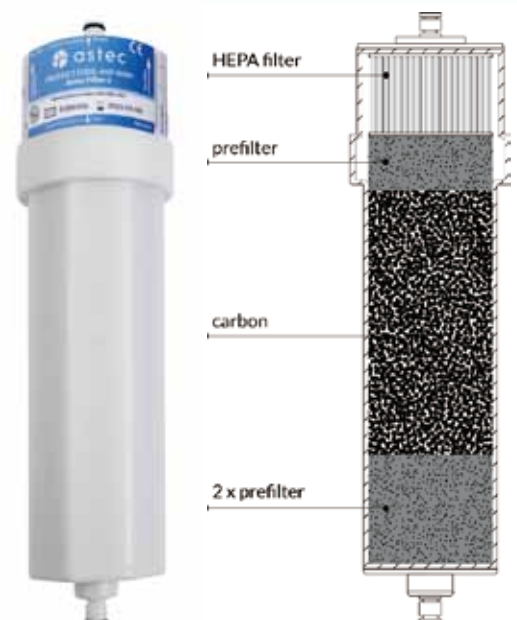

Two types of Durability

Astec Gas Line filters are available in two variants: Astec Filter 6 (six month longevity) and the longer lasting Astec Filter 4 (four month longevity).



Line Connector : AS-852085



Holder : AS-852585

High Quality

Manufactured and tested in Clean Room Environment, according to cGMP under strict rules for IVF use.

Pure Materials for IVF use

Capsule and all internal materials are safe to use in IVF laboratories.

No Chemicals Inside

Chemicals may have a negative effect on the embryo development and we have entirely eliminated the use of chemicals in our assembly and production line

Quality Control

Each filter is tested in a Clean Room Environment with an Automated Quality Testing Device, which eliminates human errors entirely.

Packing

Filters are individually packed and delivered as a set of four pieces in a sealed cardboard box which includes a certificate of conformity and instruction for use.